
CleanCare

5000 Series

Rechargeable
sonic toothbrush

PHILIPS
sonicare
the sonic toothbrush

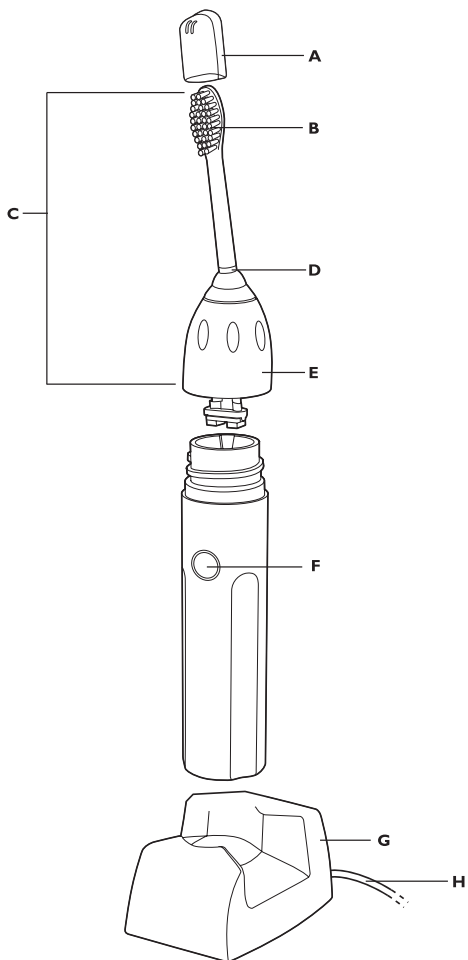ENGLISH	4
DANSK	20
ESPAÑOL	36
SUOMI	52
NORSK	67
PORTUGUÊS	83
РУССКИЙ	100
SVENSKA	117

Introduction

Congratulations on your purchase and welcome to Philips! To fully benefit from the support that Philips offers, register your product at www.philips.com/welcome.

General description (Fig. 1)

- A** Travel cap
- B** Extra-soft bristles
- C** Brush head
- D** Colour code ring
- E** Removable nut
- F** On/off button with charging light
- G** Charger with cord storage hooks
- H** Mains cord

Important

Read this user manual carefully before you use the appliance and save it for future reference.

Danger

- Keep the charger away from water. Do not place or store it over or near water contained in a bathtub, washbasin, sink etc. Do not immerse the charger in water or any other liquid. After cleaning, make sure the charger is completely dry before you connect it to the mains.

Warning

- Check if the voltage indicated on the appliance corresponds to the local mains voltage before you connect the appliance.
- Never use the charger if it is damaged in any way. If the charger, the cord or the plug is damaged, always have the charger replaced with one of the original type in order to avoid a hazard.

- If the appliance is damaged in any way (brush head, toothbrush handle, charger), stop using it. This appliance contains no user-serviceable parts. If the appliance is damaged, contact the Consumer Care Centre in your country (see chapter 'Guarantee & service').
- Do not use the charger outdoors or near heated surfaces.
- This appliance is not intended for use by persons (including children) with reduced physical, sensory or mental capabilities, or lack of experience and knowledge, unless they have been given supervision or instruction concerning use of the appliance by a person responsible for their safety.
- Children should be supervised to ensure that they do not play with the appliance.

Caution

- Do not use the toothbrush in the bath or shower.
- If you have had oral or gum surgery in the previous 2 months, consult your dentist before you use the toothbrush.
- Consult your dentist if excessive bleeding occurs after using this toothbrush or if bleeding continues to occur after 1 week of use.
- The Sonicare toothbrush complies with the safety standards for electromagnetic devices. If you have a pacemaker or other implanted device, contact your physician or the device manufacturer of the implanted device before you use the Sonicare.
- This appliance has only been designed for cleaning teeth, gums and tongue. Do not use it for any other purpose. Stop using the appliance and contact your doctor if you experience any discomfort or pain.
- If you have medical concerns, consult your doctor before you use the Sonicare.

- The Sonicare toothbrush is a personal care device and is not intended for use on multiple patients in a dental practice or institution.
- Stop using a brush head with crushed or bent bristles. Replace the brush head every 3 months or sooner if signs of wear appear.
- Do not clean the brush head in the dishwasher.
- Do not use other brush heads than the ones recommended by the manufacturer.
- Supervise children while they are brushing until they have established good brushing habits.
- If your toothpaste contains peroxide, baking soda or bicarbonate (common in whitening toothpastes), thoroughly clean the brush head with soap and water after each use. This prevents possible cracking of the plastic.

Electromagnetic fields (EMF)


This Philips appliance complies with all standards regarding electromagnetic fields (EMF). If handled properly and according to the instructions in this user manual, the appliance is safe to use based on scientific evidence available today.

Preparing for use


- 1** To store excess cord, wind it round the cord storage hooks in the bottom of the charger.

Attaching the brush head

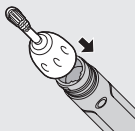- 1** To place the brush head on the handle, slide the ribs on the inside of the nut into the grooves in the top of the handle.

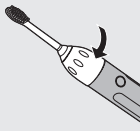
- 2** Screw the brush head nut clockwise until it is firmly tightened. If the brush head rattles during use, tighten the nut further.
- 3** Remove the travel cap from the brush head before you start brushing.

Charging

Charge the Sonicare for at least 24 hours during the first few days of use.

We advise you to keep the Sonicare in the plugged-in charger when not in use, to maintain a full battery charge.


When the battery is empty, it takes at least 24 hours to charge the Sonicare fully.

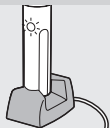- 1** Place the handle in the charger (with the on/off button facing towards the front of the charger) to charge the appliance.

Make sure the bottom of the handle sits back securely in the charger so that it makes complete contact.

- ▶ The charging light flashes slowly when the handle is in the charger to indicate that charging is in progress.
- ▶ When the Sonicare is fully charged, the charging light lights up continuously to indicate that the Sonicare is ready for use.

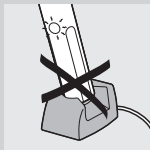Note: When the handle is placed in the charger incorrectly, the charging light may still go on, but charging does not take place.

Note: The handles of other Sonicare models do not fit in this charger and therefore will not charge effectively. Do not attempt to force a handle into this charger, as this causes damage.

Battery charge indication (specific models only)


When the battery runs low, you hear 5 beeps after the 2-minute brushing cycle and the charging light flashes rapidly. When you hear the low charge indication, there are 3 or fewer 2-minute brushings left.


Using the appliance

Brushing technique

Follow these simple steps to maximise your brushing experience.

- 1** Before you switch on the Sonicare, place the bristles of the brush head on your gumline at a 45-degree angle.

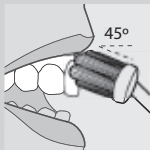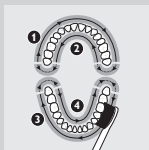
- 2** Hold the handle with a light grip and switch on the Sonicare. Gently move the bristles in a small circular motion so that the longer bristles reach between your teeth. After a few seconds, pass the bristles to the next section. Continue this routine throughout the 2-minute brushing cycle.
- 3** To maximise the effectiveness of the Sonicare toothbrush, use light pressure and let the Sonicare do the brushing for you. Do not scrub.
- 4** To make sure you brush evenly throughout the mouth, divide your mouth into 4 sections: outside upper teeth, inside upper teeth, outside lower teeth and inside lower teeth. Begin brushing the first section (outside upper teeth) and brush for 30 seconds before moving to the next section (inside upper teeth). Continue brushing the other two sections for 30 seconds each until you have completed the 2-minute brushing cycle.
- 5** After the 2-minute brushing cycle, you can briefly switch on the Sonicare again to brush the chewing surfaces of your teeth.

You can press the on/off button to switch the Sonicare on or off as desired.

Tip: To clean your mouth more completely, you can also brush your tongue with the toothbrush switched on or off, as you prefer.

Tips for specific oral care needs

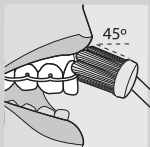- Stain removal

For even better stain removal, dedicate a few seconds of extra brushing time to areas where excess staining occurs.

- Braces

Place the brush head at a 45-degree angle between the brackets and the gumline. Using light pressure, gently move the bristles in a small circular motion so the longer bristles reach between your teeth for 2 seconds. Then gently roll the brush handle so the bristles clean the bracket and the rest of the tooth. Subsequently pass the bristles to the next tooth. Continue this procedure until you have brushed all your teeth. Please note that brush heads wear out sooner when used on braces.

- Dental restorations

Use the Sonicare on dental restorations (such as implants, caps, crowns, bridges and fillings) as you do on your natural tooth surfaces. Work the bristles round all the areas of the restoration for thorough cleaning.

- Periodontal pockets

Your dentist or dental hygienist may have indicated that you have periodontal pockets, deep pockets where plaque bacteria collect. If this is the case, spend extra time brushing these areas after you have completed the 2-minute brushing cycle to improve your gum health.

Features

Easy-start

- This Sonicare model comes with the Easy-start feature activated, which helps you get used to brushing with the Sonicare.
- The brushing power gently increases to full power over the first 12 brushings.
- The Easy-start feature is active if you hear 2 beeps immediately after you switch on the Sonicare. Once the Sonicare has reached optimal brushing power (after the 12th brushing), you no longer hear the 2 beeps at the beginning of each brushing cycle.

Note: Each of the first 12 brushings must be at least 1 minute long to move properly through the Easy-start ramp-up cycle.

To deactivate or activate the Easy-start feature:

- 1** Attach the brush head to the handle.
- 2** Place the handle in the plugged-in charger.
- 3** To deactivate the Easy-start feature, press and hold the on/off button for 5 seconds. You hear 1 beep to indicate that the Easy-start feature has been deactivated.
- 4** To activate the Easy-start feature, press and hold the on/off button for 5 seconds. You hear 2 beeps to indicate that the Easy-start feature has been activated.

Use of the Easy-start feature beyond the initial ramp-up period is not recommended and reduces


the effectiveness of the Sonicare toothbrush in removing plaque.

Make sure the Easy-start feature is deactivated and the handle fully charged when the Sonicare is used in clinical trials.

Smartimer


- The Smartimer automatically sets the brushing time to the 2 minutes that dentists and dental hygienists recommend.
- After 2 minutes, the Smartimer automatically switches off the toothbrush.
- If you want to stop brushing briefly during a 2-minute brushing cycle, press the on/off button. When you press the on/off button again, the Smartimer resumes the brushing cycle at the point where you stopped.
- The Smartimer automatically resets to the beginning of a 2-minute cycle if:
 - 1 you stop brushing for 30 seconds or longer;
 - 2 you put the handle in the charger.

Quadpacer (specific models only)

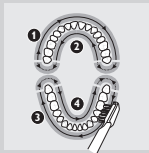- The Quadpacer 30-second interval timer ensures that you brush all sections of your mouth evenly and thoroughly. At 30, 60 and 90 seconds, you hear a short beep and notice a pause in brushing action. These signals tell you to continue with the next section of your mouth.


Cleaning

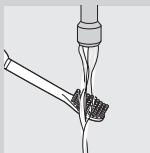Regular cleaning guarantees the best performance and a long life for the appliance.


We advise you to use a mild cleaning agent, water and a moist cloth to clean the Sonicare.

1 Rinse the brush head and bristles after every use.

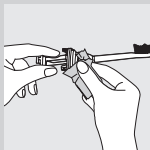
2 Clean the brush head at least once a week, including under the removable nut.

To remove the nut:

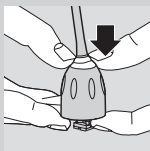
3 To remove the nut, gently press on the rubber seal with your thumbs.

Do not use sharp objects to press on the seal, as this may cause damage.

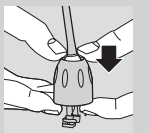4 Pull up the nut over the top of the brush head.

5 To reattach the nut, slip the nut over the brush head shaft and pull it down until it snaps home.

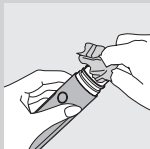6 Clean and dry the handle weekly.

7 Clean the charger periodically.

Unplug the charger before you clean it. After you have cleaned the charger, dry it completely before you put the plug back into the wall socket.

Storage

If you do not intend to use the appliance for a long time, unplug it, clean it and store it in a cool and dry place away from direct sunlight. You can store the cord by winding it round the cord storage hooks in the bottom of the charger.

Replacement

Brush head

- Replace the brush head every 3 months to achieve optimal results.
- Only use Sonicare replacement brush heads that are suitable for this Sonicare model.

Environment

- The built-in rechargeable batteries contain substances that may pollute the environment. Always remove the batteries before you discard and hand in the appliance at an official collection point. Dispose of the batteries at an official collection point for batteries. If you have trouble removing the batteries, you can also take the appliance to a Philips service centre. The staff of this centre will remove the batteries for you and will dispose of them in an environmentally safe way.


- Do not throw away the appliance with the normal household waste at the end of its life, but hand it in at an official collection point for recycling. By doing this, you help to preserve the environment.


Removing the rechargeable batteries

This process is not reversible.


To remove the batteries yourself, you need a flat-head (standard) screwdriver.

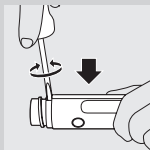
- 1** Remove the toothbrush from the charger, switch it on and let it run until it stops. Continue to switch on the Sonicare until the batteries are completely empty.

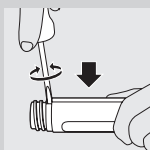
- 2** Insert a screwdriver into one of the slots that are located near the top part of the handle. Push the screwdriver blade into the handle as far as it goes and twist the blade. The housing begins to split.

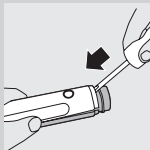- 3** Remove the screwdriver and apply the above instruction to the other slot, which is located on the opposite side of the handle. Twist the blade to split the rest of the housing.

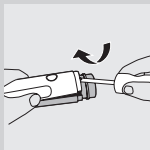- 4** Now insert the screwdriver into the middle slot, which is located on the top front part of the handle, as shown. Push the blade in as far as it goes.

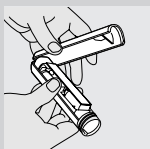
- 5** Now use the screwdriver as a lever to pry the 2 handle sections apart.

- 6** Pull the 2 handle sections completely apart.


- 7** Insert the screwdriver underneath the electronic circuit board in multiple places.

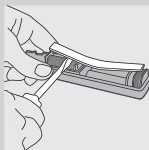- 8** Twist and pry the electronic circuit board loose until it is completely separated from the handle.

- 9** Pry the batteries loose by inserting the screwdriver on the sides of the two batteries in the locations that are shown in the figure.

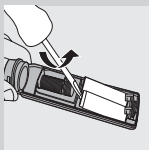- 10** Twist the screwdriver to break the glue bond round the batteries. The batteries are joined together with glue and are lifted out together.

Guarantee & service

If you need service or information or if you have a problem, please visit the Philips website at www.philips.com or contact the Philips Consumer Care Centre in your country (you find its phone number in the worldwide guarantee leaflet). If there is no Consumer Care Centre in your country, go to your local Philips dealer.

Guarantee restrictions

The terms of the international guarantee do not cover the following:

- Brush heads
- Damage caused by misuse, abuse, neglect or alterations
- Normal wear and tear, including chips, scratches, abrasions, discolouration or fading

Troubleshooting

This chapter summarises the most common problems you could encounter with the Sonicare. If you are unable to solve the problem, contact the Philips Consumer Care Centre in your country (see chapter 'Guarantee & service') or visit www.philips.com.

Problem	Solution
The Sonicare toothbrush does not work.	You may need to recharge the Sonicare. Charge the Sonicare for at least 24 hours.
	Make sure the socket is live when you recharge the appliance. The electricity supply to sockets in mirror cabinets over washbasins or shaver sockets in bathrooms may be cut off when the light is switched off.
	You may need to reattach or replace the brush head.
I feel a tickling sensation when I use the Sonicare.	You may experience a slight tickling/tingling sensation when you use the Sonicare for the first time. As you become accustomed to brushing with the Sonicare, this sensation diminishes.
The brush head rattles when I use the Sonicare.	The Sonicare should produce a gentle hum during operation. Tighten the brush head nut to eliminate excessive rattling. If the rattling does not stop, you may need to replace the brush head.
	Make sure you remove the travel cap before you switch on the appliance.
The brush head chatters against my teeth.	You can avoid this by keeping the brush head properly oriented, with the brush head placed on the gumline at an angle of 45 degrees.

Problem	Solution
Brushing seems less powerful.	You may need to recharge the Sonicare. The Easy-start feature may be active and needs to be deactivated (see chapter 'Features'). You may need to clean the brush head. You may need to replace the brush head.
The Sonicare gives off an unpleasant smell.	You may need to clean the appliance. Clean the handle, brush head and charger thoroughly (see chapter 'Cleaning'). You may need to replace the brush head.

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All rights reserved. PHILIPS and the Philips shield are registered
trademarks of KPENV. Quadpacer, Sonicare and the Sonicare logo are
trademarks of Philips Oral Healthcare, Inc. and/or KPENV.