
sonicare[®]
the sonic toothbrush

elite
e7000
e9000

PHILIPS

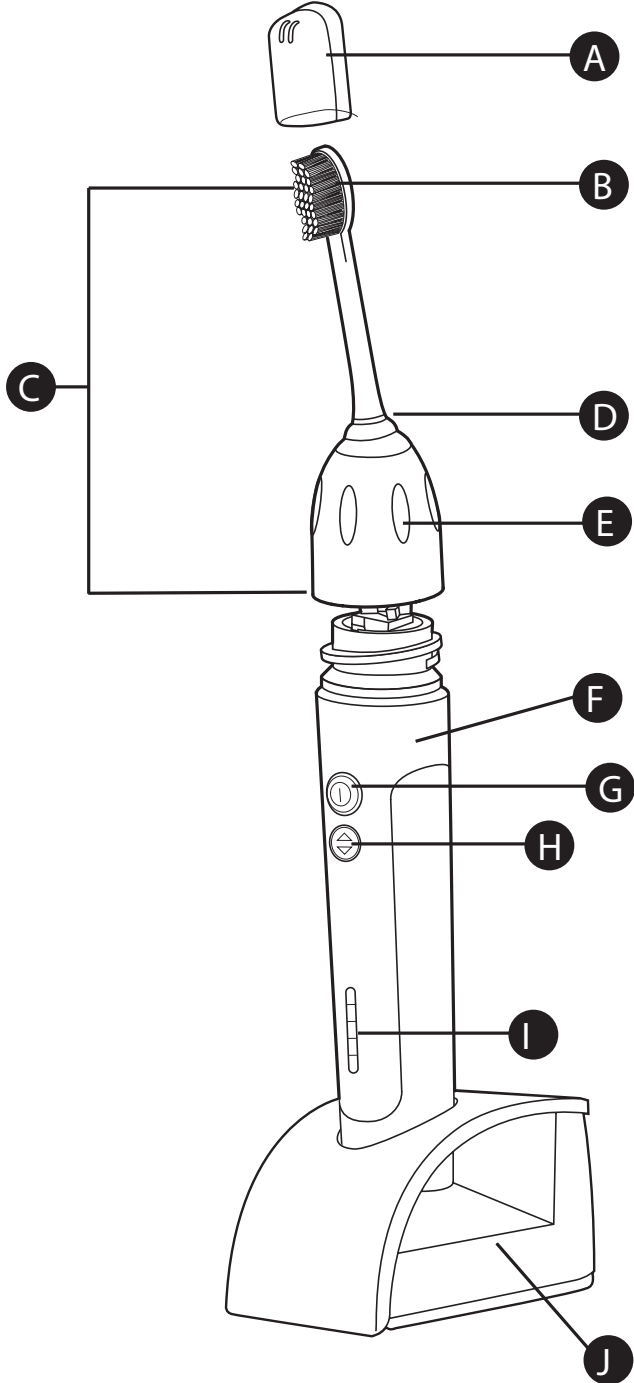Sonicare elite e7000/e9000

ENGLISH	4
DEUTSCH	15
FRANÇAIS	26
ITALIANO	37
NEDERLANDS	48

General description (Fig. 1)

- A** Travel cap
- B** Bristles
- C** Brush head
- D** Colour code ring
- E** Removable nut
- F** Handle with soft grip
- G** Power on/off button
- H** Dual Speed Control button
- I** LED charging lights
- J** Charger

Important

Read this user manual carefully before you use the appliance and save it for future reference.


- Check if the voltage indicated on the bottom or back of the charger corresponds to the local mains voltage before you connect the appliance.
- Keep the charger away from water. Do not place or store it over or near water contained in a bathtub, washbasin, sink etc. Do not immerse the charger in water or any other liquid. After cleaning, make sure the charger is completely dry before you connect it to the mains.
- Never use the charger if it is damaged in any manner. If the charger, the cord or the plug is damaged, always have the charger replaced with one of the original type in order to avoid a hazard. Do not use the charger outdoors or near heated surfaces.
- Do not use the toothbrush in the bath or shower.
- If you have had oral or gum surgery in the previous 2 months, consult your dentist before you use the toothbrush.
- Consult your dentist if excessive bleeding occurs after using this toothbrush or if bleeding continues to occur after 1 week of use.
- The Sonicare toothbrush complies with the safety standards for electromagnetic devices. If you have a pacemaker or other implanted device, consult your doctor or the manufacturer of the implanted device if you have any questions about using the Sonicare.
- This appliance has only been designed for cleaning teeth, gums and tongue. Do not use it for any other purpose. Stop using the appliance and contact your doctor if you experience any discomfort or pain.
- The Sonicare toothbrush is a personal care device and is not intended for use on multiple patients in a dental practice or institution.
- If the appliance is damaged in any way (brush head, toothbrush handle, charger), stop using it. This appliance contains no serviceable parts. If the appliance is damaged, contact the Customer Care Centre in your country (see chapter 'Guarantee & service').
- Stop using a brush head with crushed or bent bristles. Replace the brush head every 6 months or sooner if signs of wear appear.
- Do not use other brush heads than the ones recommended by the manufacturer.
- Close supervision is necessary when this product is used by or on children under the age of 10 or the disabled.

- If your toothpaste contains peroxide, baking soda or bicarbonate (common in whitening toothpastes), thoroughly clean the brush head with soap and water after each use.

Electromagnetic fields (EMF)


This Philips appliance complies with all standards regarding electromagnetic fields (EMF). If handled properly and according to the instructions in this user manual, the appliance is safe to use based on scientific evidence available today.

Preparing for use

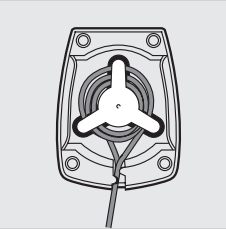
- 1 To store excess cord, wind it round the cord storage hooks in the bottom of the charger.

Attaching the brush head

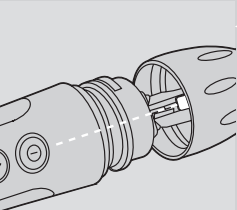
- 1 Align the front of the brush head with the front of the toothbrush handle. Match the bars on the brush head with the notches inside the top of the handle and slide the brush head onto the handle.

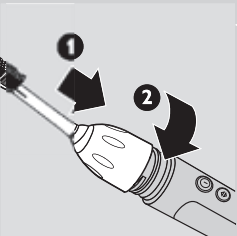- 2 Screw the brush head nut clockwise until it is firmly tightened. If the brush head rattles during use, tighten the nut further.
- 3 Remove the travel cap from the brush head before you start to brush.

Charging


Charge the Sonicare for at least 24 hours during the first few days of use. We advise you to keep the Sonicare in the plugged-in charger when not in use to maintain a full battery charge. When the battery is empty, it takes at least 24 hours to fully charge the Sonicare.

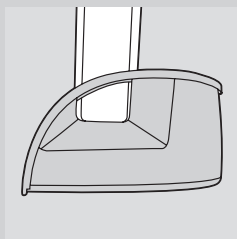1 Make sure the bottom of the handle sits back securely in the charger so that it makes complete contact.

- The handles of other Sonicare models do not fit in this charger and will therefore not charge effectively. Do not attempt to force a handle into this charger, as this causes damage.

Battery-low indicator

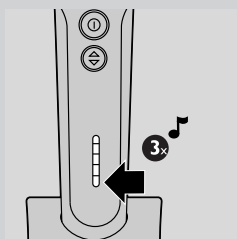When the battery runs low, you hear 3 beeps after the Smartimer has switched off the toothbrush. In addition, the bottom LED flashes yellow for 30 seconds.

When you see and hear the battery-low indication, 3 or fewer 2-minute brushings remain.

Keep the Sonicare in the plugged-in charger when not in use to maintain a full battery charge.

Multiple-LED charge indicator

When you place the handle in the charger, the green charge indicator LEDs flash in a rising pattern as the battery charges. The number of flashing LEDs indicates the level to which the battery has been charged.

When the battery is fully charged, the LEDs stop flashing and start burning continuously. About once a minute, the charge indicator LEDs go out briefly and then go on again. This indicates that the battery is maintained at full charge. The full-charge indication persists until you remove the handle from the charger.


Using the appliance

Useful tips


- The brush head comes with a convenient travel cap. Remove the travel cap from the brush head before you switch on the Sonicare.
- Rinse the brush head and bristles after each use.
- Clean the brush head weekly, including under the nut.
- A fully charged Sonicare allows you to brush your teeth for 2 minutes twice a day for approximately 2 weeks without charging. This makes it easy to travel without the charger.
- Voltage converters and plug adapters DO NOT guarantee voltage compatibility. Always check for voltage compatibility.

Brushing technique

Follow these simple steps to maximise your brushing experience.

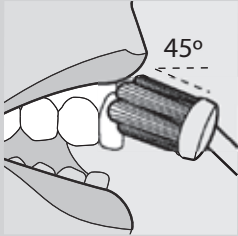- 1** Before you switch on the Sonicare, place the bristles of the brush head on your gumline at a 45-degree angle.
- 2** Hold the handle with a light grip and switch on the Sonicare. Gently move the bristles in a small circular motion so that the longer bristles reach between your teeth. After a few seconds, pass the bristles to the next section. Continue this routine throughout the 2-minute brushing cycle.
- 3** To maximise the effectiveness of the Sonicare toothbrush, use light pressure and let the Sonicare do the brushing for you. Do not scrub.
- 4** To make sure you brush evenly throughout the mouth, divide your mouth into 4 sections: outside upper teeth, inside upper teeth, outside lower teeth and inside lower teeth. Begin brushing the first section (outside upper teeth) and brush for 30 seconds before moving to the next section (inside upper teeth). Continue brushing the other two sections for 30 seconds each until you have completed the 2-minute brushing cycle.
- 5** After the 2-minute brushing cycle, you can briefly switch on the Sonicare again to brush the chewing surfaces of your teeth.

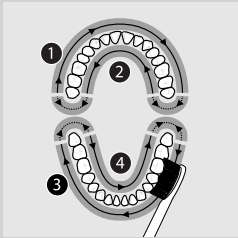- 6** Press the on/off button to switch the Sonicare on or off as needed.
 - You can also brush your tongue with the toothbrush switched on or off, as you prefer.

Tips for specific oral care needs

- Stain removal

For even better stain removal, dedicate a few seconds of extra brushing time to areas where excess staining occurs.

- Braces

Place the brush head at a 45-degree angle between the brackets and the gumline. Using light pressure, gently move the bristles in a small circular motion so the longer bristles reach between your teeth for 2 seconds.


Then gently roll the brush handle so the bristles clean the bracket and the rest of the tooth. Subsequently, pass the bristles to the next tooth. Continue this procedure until you have brushed all your teeth. Please note that brush heads wear out sooner when used on braces.

- Dental restorations

Use the Sonicare on dental restorations (such as implants, caps, crowns, bridges and fillings) as you do on your normal tooth surfaces. Work the bristles round all areas of the restoration for thorough cleaning.

- Periodontal pockets

Your dentist or dental hygienist may have indicated that you have periodontal pockets, deep pockets where plaque bacteria collect. If this is the case, spend extra time brushing these areas after you have completed the 2-minute brushing cycle to improve your gum health.

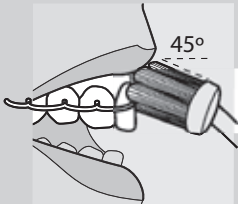Features

Easy-start

- This Sonicare model comes with the Easy-start feature activated, which helps you get used to brushing with the Sonicare.
- The brushing power gently increases to full power over the first 14 brushings.
- The Easy-start feature is active if you hear 2 beeps immediately after you switch on the Sonicare. Once the Sonicare has reached optimal brushing power (after the 14th brushing), you no longer hear the 2 beeps at the beginning of each brushing cycle.

Note: Each of the first 14 brushings must be at least 1 minute long to move properly through the Easy-start ramp-up cycle.

Note: While the Easy-start feature is activated, the Dual Speed feature (Gentle and Max speed modes) does not function. To use the Dual Speed feature immediately, deactivate the Easy-start feature.

Deactivating or activating the Easy-start feature

- 1** Attach the brush head to the handle.
- 2** Place the handle in the plugged-in charger.
- 3** To deactivate the Easy-start feature, press and hold the on/off button for 5 seconds. You hear 1 beep to indicate that the Easy-start feature has been deactivated.
- 4** To activate the Easy-start feature, press and hold the on/off button for 5 seconds. You hear 2 beeps to indicate that the Easy-start feature has been activated.


Use of the Easy-start feature beyond the initial ramp-up period is not recommended and reduces the effectiveness of the Sonicare toothbrush in removing plaque.

Make sure the Easy-start feature is deactivated and the handle fully charged when the Sonicare is used in clinical trials.

Smartimer

- This Sonicare model comes with the Smartimer feature, which automatically sets the brushing time to the 2 minutes that dentists and dental hygienists recommend.
- After 2 minutes, the Smartimer automatically switches off the toothbrush.
- If you want to stop brushing briefly during a 2-minute brushing cycle, press the on/off button. When you press the on/off button again, the Smartimer resumes the brushing cycle at the point where you stopped.
- The Smartimer automatically resets to the beginning of a 2-minute cycle if:
 - you stop brushing for 30 seconds or longer;
 - you put the handle in the charger.

Quadpacer®

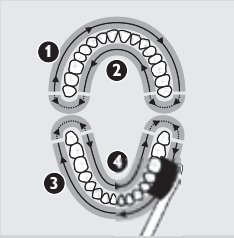- The Quadpacer 30-second interval timer ensures that you brush all sections of your mouth evenly and thoroughly. At 30, 60 and 90 seconds, you hear a short beep and notice a pause in brushing action. These signals tell you to continue with the next section of your mouth.



Programmable Quadpacer (specific types only)

The Programmable Quadpacer allows you to customise the Quadpacer in 1 of 3 ways:


- You may choose to keep the default setting of a 2-minute brushing time with the Quadpacer activated.
- You may choose to extend the brushing time by 30 seconds to 2.5 minutes with the Quadpacer activated.
- You may choose a 2-minute brushing time with the Quadpacer deactivated.

To change the Quadpacer setting:


- 1** Place the handle in the plugged-in charger.
- 2** Press the Dual Speed Control button briefly and note which LED lights up

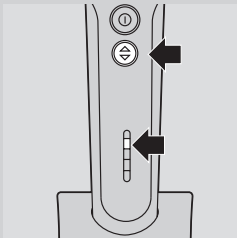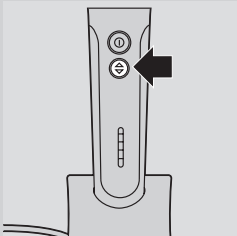- 3** Press and hold the Dual Speed Control button. Continue to hold the button while the LEDs slowly cycle through the Quadpacer settings:
 - ▶ The fourth LED burns green: the 2-minute Quadpacer is active (default setting).

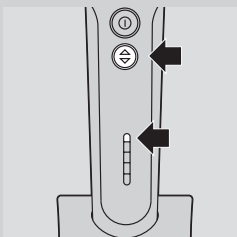- ▶ The top LED burns green: the 2.5-minute Quadpacer is active.

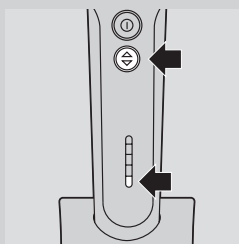► The bottom LED burns yellow: the Quadpacer is not active.

4 To select the desired setting, release the button when the LED that corresponds to your choice lights up.

Note: If you select the 2.5-minute Quadpacer, fewer brushings are available from a fully charged battery.

Dual Speed Control

The Dual Speed Control feature allows you to choose between 2 speed settings: normal speed ('Max') for optimum brushing performance or reduced speed ('Gentle') when you desire an even more gentle brushing action.

Each time you switch on the Sonicare, it starts to brush at the normal speed setting.

1 To switch from normal to reduced speed, press the Dual Speed Control button once during brushing.

2 To switch from reduced to normal speed, press the Dual Speed Control button a second time during brushing.


A distinctive two-tone beep indicates that the speed setting has been changed.

We advise you to brush at the normal speed setting to maintain optimal results.

Note: While the Easy-start feature is activated, the Dual Speed feature (Gentle and Max speed modes) does not function. To use the Dual Speed feature immediately, deactivate the Easy-start feature (see section 'Easy-start' in this chapter).


Make sure you select the normal speed setting when the Sonicare is used in clinical trials.

Cleaning

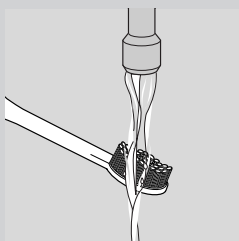
Regular cleaning guarantees the best performance and a long life for the appliance.

1 Rinse the brush head and bristles after each use. Let the brush head air-dry only.

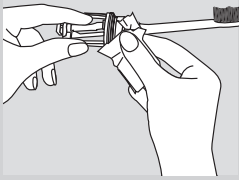
2 Clean the brush head at least once a week, including under the nut.

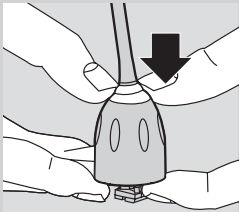
3 To remove the nut, gently press on the rubber seal with your thumbs.

Do not use sharp objects to press on the seal, as this may cause damage.

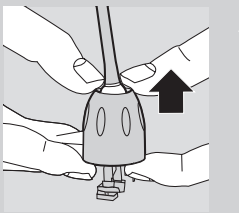- Pull the nut over the top of the brush head.

4 To reattach the nut, slip the nut over the brush head shaft and pull it down until it snaps home.

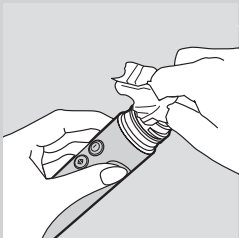5 Clean and dry the handle weekly.

6 Clean the charger periodically.

Unplug the charger before you clean it. After you have cleaned the charger, dry it completely before you put the plug back into the wall socket.

We advise you to use a mild cleaning agent, water and a moist cloth to clean the Sonicare.

Storage

1 If you are not going to use the appliance for a long time, unplug it, clean it and store it in a cool and dry place away from direct sunlight. You can store the cord by winding it round the cord storage hooks in the bottom of the charger.

Replacement

Brush head

- Replace your brush head every 6 months to achieve optimal results.
- Only use Sonicare replacement brush heads that are suitable for this Sonicare model.

Environment


- Do not throw away the appliance with the normal household waste at the end of its life, but hand it in at an official collection point for recycling. By doing this, you help to preserve the environment.
- The built-in rechargeable NiCd battery contains substances that may pollute the environment. Always remove the battery before you discard and hand in the appliance at an official collection point. Dispose of the battery at an official collection point for batteries. If you have trouble removing the battery, you can also take the appliance to a Philips service centre. The staff of this centre will remove the battery for you and will dispose of it in an environmentally safe way.

Removing the battery

This process is not reversible.


To remove the battery yourself, you need a flat-head (standard) screwdriver.

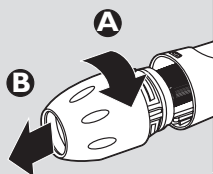
1 Remove the toothbrush from the charger, switch it on and let it run until it stops. Continue to switch on the Sonicare until the battery is completely empty.

2 Remove the nut from the brush head and place it onto the handle. Tighten the nut onto the handle (A) until the housing comes loose and then pull the nut forward (B).

The components inside the handle come out of the housing along with the nut.

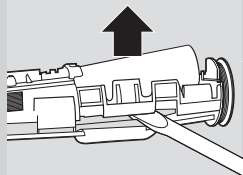3 Insert a screwdriver between the electronic circuit board and the plastic part close to each end of the battery. Twist the screwdriver to break the connection between the battery and the electronic circuit board in 2 locations.

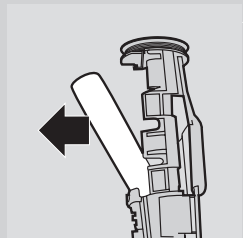4 Separate the battery from the appliance with a screwdriver.

Do not connect the appliance to the mains again after you have removed the battery.

Guarantee & service

If you need information or if you have a problem, please visit the Philips website at www.philips.com or contact the Philips Customer Care Centre in your country (you find its phone number in the worldwide guarantee leaflet). If there is no Customer Care Centre in your country, go to your local Philips dealer or contact the Service Department of Philips Domestic Appliances and Personal Care BV.

Guarantee restrictions

The terms of the international guarantee do not cover the following:

- Brush heads
- Damage caused by misuse, abuse, neglect or alterations.
- Normal wear and tear, including chips, scratches, abrasions, discolouration or fading.

Frequently asked questions

This chapter summarises the most common problems you could encounter with the Sonicare. If you are unable to solve the problem, contact the Philips Customer Care Centre in your country (see chapter 'Guarantee & service') or visit www.philips.com.

Question

Answer

Why doesn't the Sonicare toothbrush work?

You may need to recharge the Sonicare. Charge the Sonicare for at least 24 hours.

Make sure the socket is live when you recharge the appliance. The electricity supply to sockets in mirror cabinets over washbasins or shaver sockets in bathrooms may be cut off when the light is switched off.

Why doesn't the brush head operate while the motor runs?

You may need to reattach or replace the brush head.

Why do I feel a tickling sensation when I use the Sonicare?

You may experience a slight tickling/tingling sensation when you use the Sonicare for the first time. As you become accustomed to brushing with the Sonicare, this sensation diminishes.

Why does the brush head rattle when I use the Sonicare?

The Sonicare should produce a gentle hum during operation. Tighten the brush head nut to eliminate excessive rattling. If the rattling does not stop, you may need to replace the brush head.

Make sure you remove the travel cap before you switch on the appliance.

The brush head chatters against my teeth. What can I do about this?

To avoid this, keep the brush head properly oriented, with the brush head placed on the gumline at an angle of 45 degrees.

Why does the toothbrush seem less powerful?

You may need to recharge the Sonicare.

The Easy-start feature may be active. Deactivate the Easy-start feature (see chapter 'Features').

Question	Answer
	You may need to clean the brush head.
	You may need to replace the brush head.
	The Dual Speed Control button may be set to reduced speed. Set the button to normal speed. You may need to clean the appliance. Clean the handle, brush head and charger thoroughly (see chapter 'Cleaning').
The Sonicare gives off an unpleasant smell. What can I do about this?	You may need to replace the brush head.
	You may need to clean the appliance. Clean the handle, the brush head and the charger thoroughly (see chapter 'Cleaning').
The yellow LED flashes. What does this mean?	The flashing yellow LED indicates that fewer than 3 brushings remain. Recharge the Sonicare for at least 24 hours.
The appliance has an extra 30 seconds of brushing time. What is this for?	The Programmable Quadpacer has been set to a 2.5-minute brushing cycle instead of the default 2-minute brushing cycle. This setting provides additional time for brushing troublesome areas where plaque collects and gum disease starts. We advise you to use the extra 30 seconds of brushing time to rebrush troublesome areas, massage your gums or brush your tongue (see chapter 'Features').

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Philips Oral Healthcare, Inc.
35301 SE Center Street
Snoqualmie WA 98065
U.S.A.

Outside North America:
+1-425-396-2000
www.philips.com

Printed in U.S.A.
Printed on 100% recycled paper

Protected by U.S. and international patents. Other patents pending.

© 2006 Philips Oral Healthcare, Inc. All rights reserved. PHILIPS and the Philips shield are registered trademarks of Koninklijke Philips Electronics N.V. Quadpacer, Sonicare, Elite and the Sonicare logo are trademarks of Philips Oral Healthcare, Inc. and/or KPEENV.